

100% Ethylene Oxide sterilization



The Office Sterilization Solution

The most popular low temperature sterilization system keeps getting better! Using less than 18 grams of Ethylene Oxide per cycle, Anprolene® gently sterilizes your delicate items that heat and moisture can damage.

Andersen's unique bag and ampoule method allows for load size flexibility.

Easy installation requires only a dedicated 110V outlet and exhaust line to the outside.

► Unmatched Sterilization

Flexibility: Anprolene's 100% EtO System sterilizes at room temperature without external vacuum lines or steam generator. Protect your most delicate instruments, scopes, drills, fiber optics, cameras, and plastic or rubber items... Literally everything but food, drugs or liquids can be sterilized with Anprolene.

► **Sterilization and Aeration in the same cabinet:** The Anprolene system utilizes an advanced purge system that aerates products for two hours at the end of the cycle. Active aeration continues indefinitely until products are removed.

► **Unique Load Size Flexibility:** The Anprolene system uses unique plastic sterilization liner bags and individual gas ampoules. We offer two different ampoule/bag combinations to match the size of your sterilization load. You never use more gas than you have to!

► **Proven Reliability:** Anprolene has been the choice of healthcare professionals for more than 30 years. Anprolene is used in over 5,000 facilities worldwide.



AN71 Refill Kit



AN93 AirScan



AN87 Dosimeter

STARTER KIT INCLUDES

| Units | Item | Reorder | Quantity |
|----------|---|---------|----------|
| 25 each | 5 cc ampoules and sterilization liner bags | AN71 | 1 case |
| 20 each | 20 cc ampoules and sterilization liner bags | AN79 | 1 case |
| 50 each | Humidichip relative humidity stabilizer chips | AN1071 | 2 jars |
| 1 each | Humiditube holder | AN1072 | 1 each |
| 50 each | Dosimeter chemical indicators | AN87 | 2 boxes |
| 400 each | EtO exposure indicator strips - chemical indicators | AN85 | 2 boxes |
| 22 each | Steritest biological indicators | AN80 | 2 boxes |
| 1 each | AirScan test kit (STEL and TWA) | AN93 | 1 kit |
| 1 each | Installation Kit | AN74.01 | 1 each |

* Sterilizer must be purchased with a starter kit.



sterilization redefined

www.anpro.com
Andersen Products, Inc.
Health Science Park • 3202 Caroline Drive
Haw River, NC 27258-9564 USA
800-523-1276 • fax 336-376-8153

www.andcal.com
H. W. Andersen Products of California, Inc.
Health Science Park • 3151 Caroline Drive
Haw River, NC 27258-9789 USA
800-524-3455 • fax 336-376-3088

Distributed in Europe by:
H.W Andersen Products, Ltd.
Davy Road • Clacton-on-Sea
Essex CO15 4XA UK
1255-428-328


Manufactured by:
Andersen Sterilizers, Inc.
3154 Caroline Drive
Haw River, NC 27258 USA

Sterilizer technical specifications

- ▶ **Anprolene uses less Gas:** The Anprolene system uses less than 18 grams of EtO per cycle. The small quantity of gas used means that the Anprolene system generally meets State and U.S.E.P.A. emission requirements without the need for additional abatement equipment.
- ▶ **Sterility Assurance:** The Anprolene System uses a standard range of sterility assurance indicators to confirm the success of the sterilization cycle. Steritest® is a Biological Indicator (BI), the classic test of sterility. BI's require a minimum 48 hr. culture period. Immediate results are provided by Dosimeter®, a Chemical Indicator, and a rigorous test for every cycle. Anprolene Exposure Indicators are used for each item within a "sterilization load" to ensure individual exposure.
- ▶ **Ease of Installation:** New high pressure, low volume pumps simplify installation, requiring only a 7/8" reinforced PVC exhaust hose to the outside. No water or vacuum line required! Wall mount installation kit included with the sterilizer.

Specifications: □ AN74i

Model AN74.67
Interior dimensions: 19.5"(L) x 13.75"(W) x 13.75"(H)
Exterior dimensions: 21"(L) x 14"(W) x 20.5"(H)
Total weight: 50.5 lbs.
Power Supply: 115v +/- 10%, 60 Hz/150 Watts
Installation Category II
Manufactured by: Andersen Sterilizers, Inc.
Distributed by: Andersen Products, Inc.

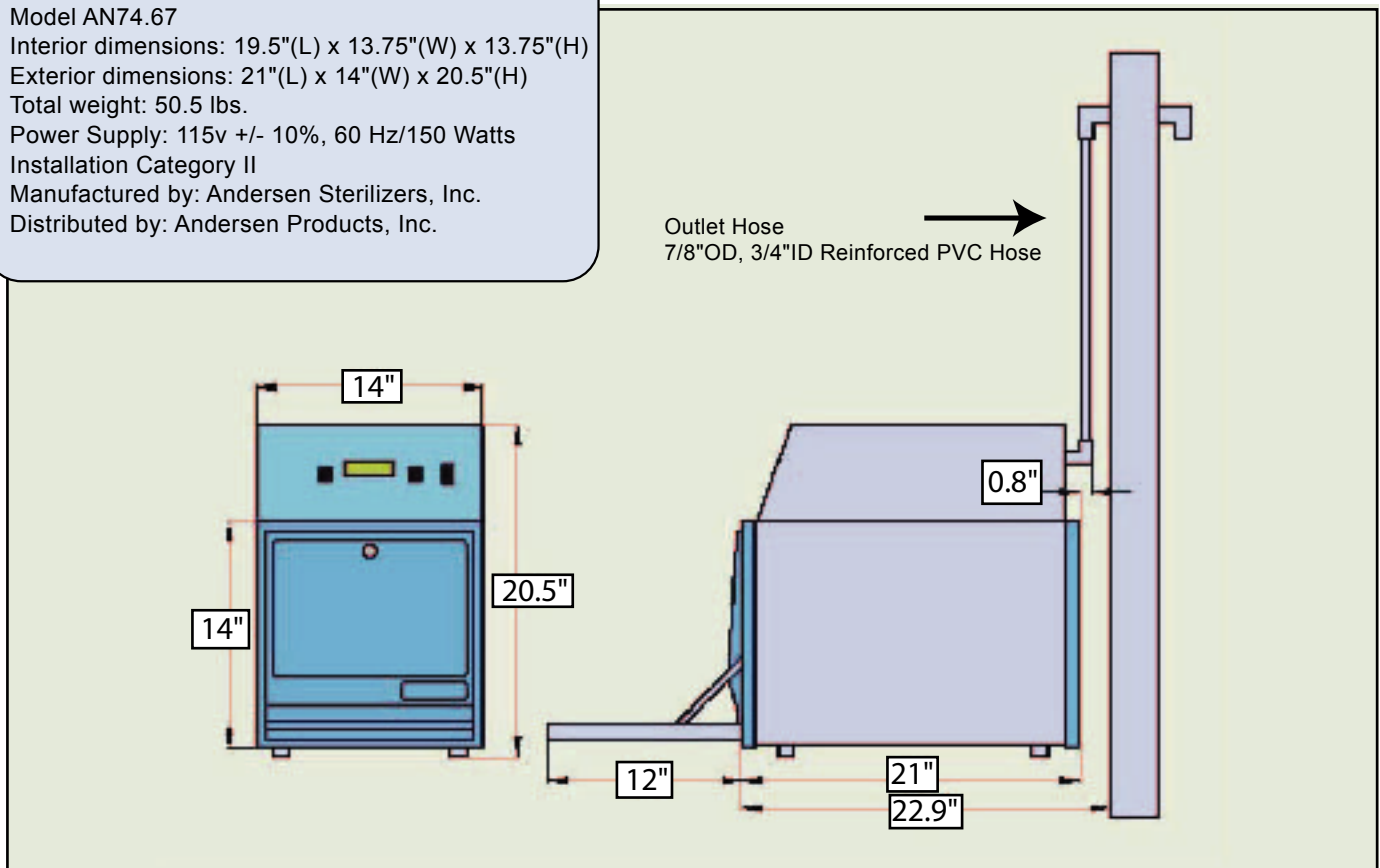

For more information about Anprolene® Sterilization Kits and prices call customer service at 800-523-1276, or visit www.anpro.com

