

# 80% less gas than the competition.

*With Andersen's match load to cartridge method you never have to worry about tanks of EtO. Andersen Products offers cutting-edge technologies and over 30 years of expertise in sterilization with the only system that offers multiload capability.*

► **Match Load Volume:**

Two bag sizes available.

The AN1005 is 18" x 24" bag size with only 5 grams of EtO used.

The AN1006 is 22" x 36" bag size with only 11 grams of EtO used.

► **Cost Effective Kits include:**

25 EOGas 100% EtO cartridges

25 Humidichip® relative humidity stabilizing chips


25 EOGas™ Dosimeters®

25 corresponding size Sterilization Bags

► **Design:**

Its robust design and polyurethane body make the EOGas™ cartridge one of the most reliable vehicles for releasing the ethylene oxide gas.

Operator safety is insured by the double protection of the trigger via a plastic lock and a safety seal.

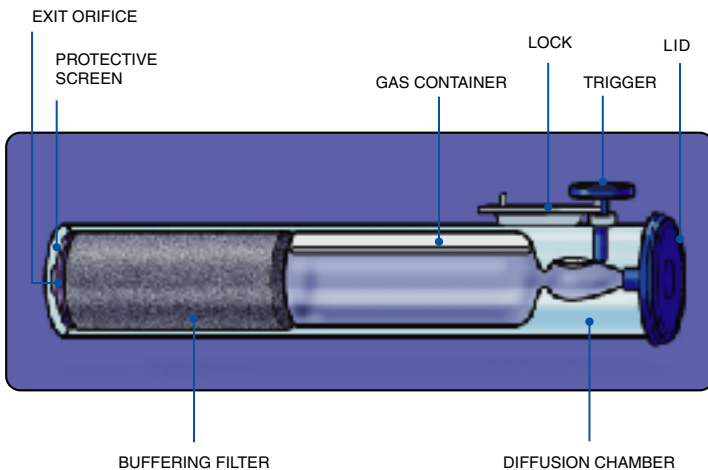

*Two sizes of sterilant cartridges and corresponding bags simplify sterilization. Process loads of various sizes simultaneously. Avoid restrictions of single chamber batch processing.*



*The gas goes through 3 barriers before it reaches the venting mechanism.*



*Double protection trigger.*

092603